

# XL-700™


## ( TALKING SCALE )



ENGLISH



**Capacity**  
320kg x 0.1kg  
700lb x 0.2lb

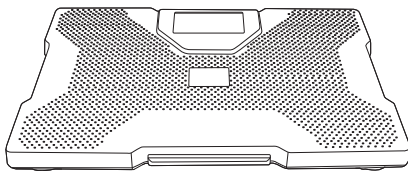

Thank you for purchasing the MyWeigh®XL700™ talking bathroom scale. Please read all operating instructions carefully before use. This electronic scale is a precision instrument. With normal care and proper treatment, it will provide years of reliable service. For more information please visit [www.myweigh.com](http://www.myweigh.com).

### OPERATION INSTRUCTIONS

#### Inserting batteries

Insert 4 xAA batteries into the battery compartment on the bottom (underside) of the scale.

#### Selecting language

After showing "0.0", press and hold the [kg/lb/st] button for 7 seconds. You will hear "Hello" of the existing language. Press the [kg/lb/st] button again to select the language. You will hear "Hello" in the other languages. To confirm the selection, press and hold the [kg/lb/st] button until it says "Goodbye".

#### Weighing

- Tap the Activation button at the front of the scale with your foot. The scale will say "Hello".
- Once the scale is stable, it'll say "It's ready" and show "0.0 kg" (or 0.0 lb, 0.0 st). Then you may step onto the measurement platform again. Your weight will be displayed and read out once the reading is steady.
- You may step off the platform and the reading will hold for 10 secs. After that it'll say "Good-bye" and turn off automatically.

#### Power off


The scale will shut down automatically after no operation is detected for about 10 seconds.

#### Switching on/off the talking function

Use the switch on the bottom of the platform to switch on/off the speaker function.

#### Changing the Unit

Turn the Scale ON, "0.0" will show on the display", use the [kg/lb/st] button under the scale to toggle the desired unit (lb, kg or st).

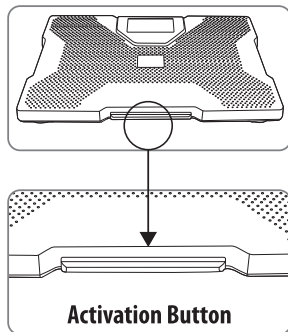

### IMPORTANT NOTES

\*\*\*If the Scale ever malfunctions please try installing new batteries, sometimes low batteries or defective batteries may cause inaccurate readings or other errors  
- When the battery runs very low, the warning icon will appear on the display. When this happens, be sure to replace the batteries immediately on the underside of the scale.

- Use the unit on hard, flat, dry and secure surface only.
- Place the scale in an area free from high temperature or high humidity. To clean, use a dampened cloth and mild detergent.
- Avoid excessive impacts or vibrations to the unit.
- If you do not intend to use this unit for a long period of time, remove the batteries before storage.
- Do not attempt to disassemble the unit, otherwise the warranty will become void.
- The unit is intended for home and personal use only (not legal for commerce).
- This scale must only be used on a hard surface such as tile, wood linoleum or short carpet.
- This and most bathroom scales will not perform properly on thick-carpeted surfaces.

<b>Capacity</b>	320kg x 0.1kg 700lb x 0.2lb	<b>Units</b>	kg, lbs, st
<b>Deviation</b>	+/- 1%		
<b>Temperature range</b>	Optimum 15-20°C		
<b>Power source</b>	4 x AA		

